Management of Diabetic Peripheral Neuropathy & PVD*

Nurite® MEDIC SHOES

Med Device for Happy Feet

- US FDA Approved & CE certified
- Proven by Clinical Trials
- Based on the Reflexology Principles
- Improves Peripheral Blood Circulation
- Provides relief from Neuropathy and PVD symptoms e.g. Pain in feet, Pain wakes at Night, Tingling, Burning

Imported from

Medic Shoes Ltd., Israel

Imported By

Health 'N' U Therapeutics Pvt. Ltd.

Room No. 7, Basement B-1,
Centaur House, CTS.No.4090,
Shanti Nagar Ind. Estate, Vakola,
Santacruz (E), Mumbai-400 055,
Maharashtra.

For More information log on to
<http://www.medicshoes.com>

Marketed By

Centaur Pharmaceuticals Pvt. Ltd.
Centaur House,
CTS. No. 4090, Shanti Nagar
Ind. Estate, Vakola, Santacruz
(E), Mumbai- 400 055.

Introducing

1st
Time in
India

US FDA Approved

Medical Device for Happy Feet

Saksham
Being Able

Centaur Pharmaceuticals Pvt. Ltd.

Centaur House, Near Grand Hyatt, Shantinagar, Vakola, Santacruz-East, Mumbai- 400 055 INDIA
Tel: 022-6649 9100 Fax : 022-6649 9108 / 112 Email : centaur@centaurlab.com

Nurite® MEDIC SHOES

Provide immediate relief from symptoms of pain caused by nerve damage in Diabetes and Disease of Blood vessels in Lower Parts of Leg (Peripheral Vascular Disease -PVD)

Benefits

- Easy to use and safe Medical Device.
- Relieves pain, Awakening in night due to pain, Tingling and Numbness.
- Reduces burning sensations.


About Nurite® MEDIC SHOES

- Nurite Medic Shoes is Medical Device which improves peripheral blood supply in feet.
- Nurite Medic Shoes were developed according to Reflexology Principles.
- Nurite Medic Shoes has raised areas within to provide pressure on specific points in your feet.
- Ergonomic straps will help you tighten the shoes & make reflexology treatment possible while sitting or lying down.
- Nurite Medic Shoes are operated by battery powered remote control.
- Nurite Medic Shoes treatment is safe & Provides medical device addition to your current medical therapy.

Med Device for Happy Feet

Nurite® MEDIC SHOES User Manual

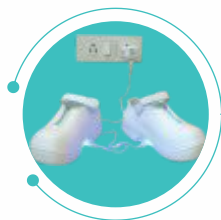Step 1:
Charging for 4 Hrs

- Connect charger as shown in image.
- Charge for 4 Hours (Blue light will turn green after complete charging)
- After green light appears, Nurite Medic Shoes are ready for use.

Step 2:
Start by pressing button

- Press the side button to start
- When Red light turns green, then Nurite Medic Shoes are ready to operate.

Step 3:
Tightening the strap

- Place both your feet inside the shoes.
- Tighten the ergonomic strap to fit properly.

Step 4:
Can be operated in Sitting
or Lying down position


- You can operate the shoes while sitting or lying down

Med Device for Happy Feet

Step 5:

Operated with Remote control

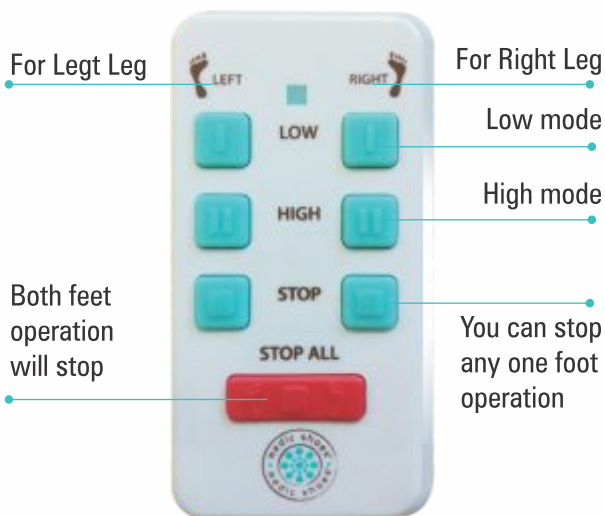- Press right side button to operate right foot shoe & Press left side buttons to operate left foot shoe.
- You can start therapy with low mode button.
- You can increase speed of therapy with high mode button.
- You can stop any of the Shoe by pressing stop button on left or right side.
- You can stop both shoes at a time by pressing red button.

Step 6:

Reflexology therapy should be taken for period of 15-30 minutes/day.

Step 7:

You should continue therapy for Minimum 30 days. For long term relief, you should continue therapy beyond 30 days.

Step 8:

Preferable time for Therapy- Before Bedtime.

Step 9:

Nurite Medic Shoes are available in 4 different sizes for Men & Women .

MEN	WOMEN
4	6
6	8
9	11
11	13

- You can complete 1 cycle for 15 minutes
- Nurite Medic Shoes will automatically stop after 15 minutes
- To reactivate Nurite Medic Shoes, You should re-press the button for each shoe.
- After three consecutive activations, Nurite Medic Shoes will automatically enter in hibernation mode and red light will blink.
- Hibernation mode will last for 1 Hour.

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